

ACCESS



# HomeGlide

The HomeGlide stairlift delivers a solution for all straight staircase requirements and budgets. Designed for effortless simplicity with superior comfort and value.



## Slim Rail

The aluminium rail has a slim design giving a modern and clean look. The HomeGlide is fitted onto your staircase and not your wall. This allows the lift to be installed quickly, efficiently and with minimal disruption to your home.


## Seatbelt

A retractable seatbelt keeps you secure while using your HomeGlide stairlift. The integrated sensor ensures the seatbelt is fastened securely before travel, and the stainless steel buckle is designed to be operated with just one hand.

## Discreet design

Engineered with modern manufacturing techniques and technology, the HomeGlide stairlift range can be fitted to almost all straight staircases as narrow as 740mm (29").

**MED**m&bility  
[www.albertastairlifts.com](http://www.albertastairlifts.com)

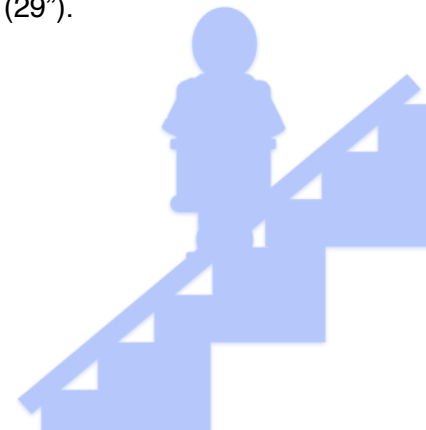

# Technical Specification

Certified to the highest standards



|                             |   |
|-----------------------------|---|
| <b>Driving Speed</b>        | 0.12 m/s (23.6 ft/min)  |
| <b>Drive Type</b>           | Rack and pinion   |
| <b>Motor Wattage</b>        | 280 W   |
| <b>Weight Capacity</b>      | 160kg (25st) at max 45°<br>136kg (21st) at max 53°<br>145kg (22st) with hinge |
| <b>Batteries</b>            | 2 × 12V, 7Ah (indoor)   |
| <b>Operation</b>            | Joystick  |
| <b>Staircase Incline</b>    | 28° up to 53° (up to 45°<br>for use in 160kg range)                           |
| <b>Certification</b>        | Machine Directive<br>2006/42/EC EN 81-<br>40:2020 UKCA                        |
| <b>Model Code</b>           | RP00-ST   |
| <b>Automatic Stop</b>       | Yes   |
| <b>Footrest</b>             | Yes   |
| <b>Retractable Seatbelt</b> | Yes   |

Elevate your home's accessibility with Access BDD's innovative HomeGlide stairlift, designed for effortless, seamless movement and superior comfort.



## Folding Rail (option)

If you have an obstruction at the foot of the staircase, you can choose to have an automatically powered hinge installed. Once you have used your stairlift the lower section of the rail will fold away, leaving the bottom of your staircase clear.



## Joystick

With a comfortable soft-touch grip, the integrated joystick is user-friendly and requires minimal effort to take you up and down your staircase with ease.



## Call Device

Our optional call device ensures you stay connected with friends and family, with up to three different telephone numbers stored and contactable directly from your stairlift.

**MED**m&bility

MEDmobility Calgary  
Unit 1, 6025 12 Street S.E.  
Calgary, AB T2H 2K1

Phone: 403-984-6676

MEDmobility Head Office  
5711 - 104 St. NW  
Edmonton, Alberta T6H 2K1

Phone: 780-437-3300